



OmniTaper® Implant System

Surgical product catalog

All Dentsply Sirona Instructions for Use are available on www.dentsplysirona.com/ifu

All products may not be regulatory cleared/released/licensed in all markets. Please contact the local Dentsply Sirona sales office for current product assortment and availability.

To improve readability for our customers, Dentsply Sirona does not use ® or ™ in body copy. However, Dentsply Sirona does not waive any right to the trademark and nothing herein shall be interpreted to the contrary.

Product illustrations are not to scale.

CONTENTS

1. Implant packaging	4
Color coding	4
2. Implants	5
3. Cover Screw	6
4. TempBase Cap	6
5. Healing abutments	6
6. Instruments	9
Initial drills	9
Drills	9
Taps	10
Implant drivers	11
Planning aids	12
Torque Wrench	14
Driver	14
Surgical Trays	15
Washtray & Overlay	17
Sterilization Container & Reusable Filter	18
7. Guided Surgery Instruments	19
Punch	19
Drills	19
Taps	21
Sleeves	22
Driver	22
Sleeves for Driver	22
Stabilization Abutment	23
Positioning Aid	23
Surgical Trays	23
8. Material specifications	26
9. Index	27

1. Implant packaging

OmniTaper implants are supplied in a double-blister package with an outer carton.



Outer box

- Specific color coding, diameter and prominent implant length information on side labels
- Data matrix code on two sides
- Stackable, with all important product information visible
- Instructions For Use (IFU) available electronically (eIFU): ifu.dentsplysirona.com



Outer blister

- Peel-off label with batch code for convenient treatment documentation



Inner blister

- Inner sterile package
- Implant shuttle with implant and cover screw
- Peel-off label with batch code for convenient treatment documentation

Color coding

The OmniTaper implant is available in five diameters and six lengths. The color coding makes it easy to identify the correct connection and to select the right prosthetic components.

Implant Ø	3.0mm	3.4mm	3.8mm	4.5mm	5.5mm
Connections	XS	S	M	L	XL
Lengths	–	–	8mm	8mm	8mm
	–	9.5mm	9.5mm	9.5mm	9.5mm
	11mm	11mm	11mm	11mm	11mm
	13mm	13mm	13mm	13mm	13mm
	15mm	15mm	15mm	15mm	15mm
	–	18mm	18mm	18mm	–

2. Implants

XS OmniTaper EV Ø3.0



Length mm	11	13	15
Order No.	68011140	68011141	68011142

S OmniTaper EV Ø3.4




Length mm	9.5	11	13	15	18
Order No.	68011143	68011144	68011145	68011146	68011147

M OmniTaper EV Ø3.8



Length mm	8	9.5	11	13	15	18
Order No.	68011148	68011149	68011150	68011151	68011152	68011153

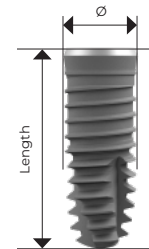

L OmniTaper EV Ø4.5



Length mm	8	9.5	11	13	15	18
Order No.	68011154	68011155	68011156	68011157	68011158	68011159

XL OmniTaper EV Ø5.5



Length mm	8	9.5	11	13	15
Order No.	68011160	68011161	68011162	68011163	68011164

OmniTaper EV Implant

Titanium¹, Titanium alloy², sterile

The OmniTaper implant is designed for replacing one or more missing teeth in order to restore chewing function.

- Endosseous section with OsseoSpeed surface
- Compatible with EV Prosthetics assortment
- Restoration: Various possibilities for type of fixation e.g., screw-, cement-, friction- and attachment-retained solutions. For one-position-only, six-position and non-indexed abutments
- Clinical use: one-stage or two-stage surgical procedures
- Premounted multifunctional TempBase that serves as a placement head, index impression and basis for temporary restorations
- Color coded
- Contains: Implant, Placement head (TempBase), Cover Screw
- Instruments needed: Implant Driver TempBase or optional Implant Driver EV

3. Cover Screw

Cover Screw EV

Titanium alloy², sterile

The Cover Screw is used for sealing the implant interface.

- Clinical use: two-stage surgical procedures
- Color coded

Cover Screw EV

					
	XS	S	M	L	XL
Order No.	68010001	68010002	68010005	68010006	68010007

4. TempBase Cap

TempBase Cap

Polypropylene²¹, non-sterile

The TempBase Cap is used for index registration during surgery and basis for temporary restorations.

TempBase Cap

					
	XS	S	M	L	XL
Order No.	68010013			68010014	68010015


5. Healing abutments

HealDesign™ EV

Titanium alloy², sterile

The HealDesign EV is used for soft tissue sculpturing during the healing phase.

- Marked with diameter and height
- Clinical use: one-stage or two-stage surgical procedure
- Round shapes are indicated for all positions in the mouth
- Triangular shapes for anterior implant sites to mimic the incisors and canines
- Triangular HealDesign EV is a two-piece abutment
- Color coded

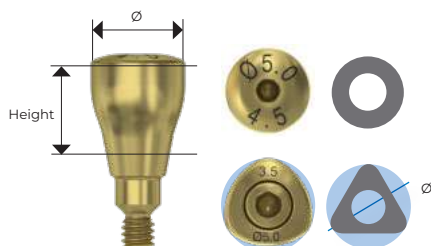

XS HealDesign EV

				
	Ø	Ø	Ø	Ø
Ø mm	3.5	3.5	3.5	3.5
Height mm	2.5	3.5	4.5	6.5
Order No.	68013001	68013002	68013003	68013004

S HealDesign EV


Ø mm	4	4	4	4	5	5	5	5	6	6	5
Height mm	2.5	3.5	4.5	6.5	2.5	3.5	4.5	6.5	3.5	4.5	3.5
Order No.	68013005	68013006	68013007	68013008	68013009	68013010	68013011	68013012	68013013	68013014	68013039

M HealDesign EV


Ø mm	5	5	5	5	6.5	6.5	6.5	6.5	5	6.5
Height mm	2.5	3.5	4.5	6.5	2.5	3.5	4.5	6.5	3.5	3.5
Order No.	68013015	68013016	68013017	68013018	68013019	68013020	68013021	68013022	68013040	68013041

L HealDesign EV


Ø mm	5	5	5	5	6.5	6.5	6.5	6.5	7.5	7.5	5	6.5
Height mm	2.5	3.5	4.5	6.5	2.5	3.5	4.5	6.5	3.5	4.5	3.5	3.5
Order No.	68013023	68013024	68013025	68013026	68013027	68013028	68013029	68013030	68013031	68013032	68013042	68013043

XL HealDesign EV



Ø mm	6.5	6.5	6.5	6.5	7.5	7.5
Height mm	2.5	3.5	4.5	6.5	3.5	4.5
Order No.	68013033	68013034	68013035	68013036	68013037	68013038

Healing Abutment EV

Titanium alloy², sterile

The Healing Abutment EV is used for soft tissue sculpturing.

- Laser-etched bands for gauging
- Marked with diameter
- Color coded

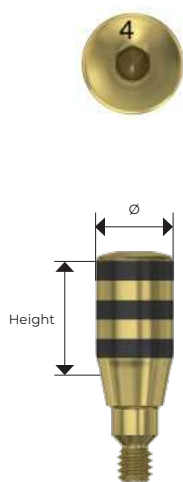

XS Healing Abutment EV

Ø mm	4	4	4	4
Height mm	2	3	4	6
Order No.	68013044	68013045	68013046	68013047

S Healing Abutment EV

Ø mm	4	4	4	4
Height mm	2	3	4	6
Order No.	68013048	68013049	68013050	68013051

M Healing Abutment EV

Ø mm	4	4	4	4
Height mm	2	3	4	6
Order No.	68013052	68013053	68013054	68013055

L Healing Abutment EV

Ø mm	4.3	4.3	4.3	4.3
Height mm	2	3	4	6
Order No.	68013056	68013057	68013058	68013059

XL Healing Abutment EV

Ø mm	4.8	4.8	4.8	4.8
Height mm	2	3	4	6
Order No.	68013060	68013061	68013062	68013063

6. Instruments

Guide Drill



Precision Drill



Initial drills

Guide Drill

Stainless steel³, sterile

The Guide Drill is used for marking and to create a starting point.

- Multiple use

Precision Drill

Stainless steel³, sterile

The Precision Drill is used for marking and to create a starting point.

- Single use

All drills, except the Precision Drill, can be used for maximum ten cases. They should be carefully cleaned and sterilized after each surgery and replaced as soon as a decrease in their cutting efficiency is observed.

The Precision Drill is an extremely sharp drill and should be handled with caution to avoid injury.

Drills

OmniTaper Drill

Stainless steel³, TiN-coated⁴, sterile

The OmniTaper Drill is used for implant site preparation.

- Two-edged profile
- Multiple use
- Depth marked
- Color coded

OmniTaper Drill short




Ø mm	2.0	3.0	3.4	3.8	4.5	5.5
Color coding						
Order No.	68015098	68015099	68015100	68015101	68015102	68015103

OmniTaper Drill long



Ø mm	2.0	3.0	3.4	3.8	4.5	5.5
Color coding						
Order No.	68015086	68015087	68015088	68015089	68015090	68015091

Drill short 8-15mm
Drill long 8-18mm

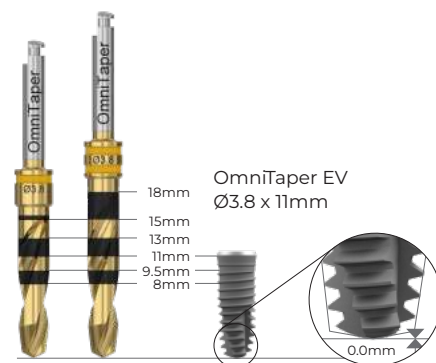

The effective drill depth matches the implant length.

OmniTaper Single Patient Drill

Stainless steel³, ABS⁹, sterile

The OmniTaper Single Patient Drill is used for implant site preparation.

- Single use
- Depth marked
- Color coded
- Available in short L 8-15mm

OmniTaper Single Patient Drill




Ø mm	2.0	3.0	3.4	3.8	4.5	5.5
Color coding						
Order No.	68015080	68015081	68015082	68015083	68015084	68015085

OmniTaper Crestal Drill

Stainless steel³, TiN-coated⁴, sterile

The OmniTaper Crestal Drill is used for crestal implant site preparation according to the bone quality once the planned implant diameter has been reached.

- Multiple use
- Depth marked
- Color coded

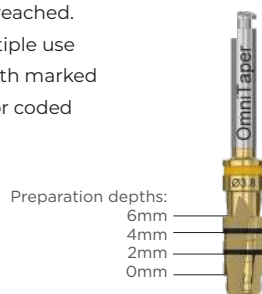

OmniTaper Crestal Drill



Ø mm	3.0	3.4	3.8	4.5	5.5
Color coding					
Order No.	68015092	68015093	68015094	68015095	68015096


Taps

OmniTaper Tap

Stainless steel³, sterile

The OmniTaper Tap is used for additional tapping in cortical bone (D I).

- Multiple use
- Color coded
- Instruments needed: contra-angle with HXSS connection or torque wrench with surgical driver handle

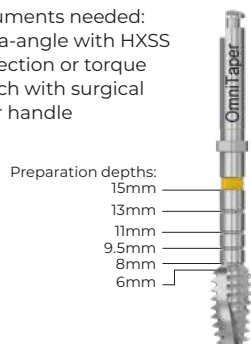

OmniTaper Tap



Ø mm	3.0	3.4	3.8	4.5	5.5
Color coding					
Order No.	68015104	68015105	68015106	68015107	68015108

Implant Driver TempBase



Order No.	68015129	68015130
-----------	----------	----------

Implant Driver TempBase O-Ring



Order No.	68015131
-----------	----------

Implant Driver EV short



Ø mm	3.4	3.8	4.5	5.5
Color coding				
Order No.	68015188	68015190	68015192	68015194

Implant Driver EV long



Ø mm	3.0	3.4	3.8	4.5	5.5
Color coding					
Order No.	68015187	68015189	68015191	68015193	68015195

Implant drivers

Implant Driver TempBase

Stainless steel¹³, EPDM¹³, non-sterile

The Implant Driver TempBase is used for installing OmniTaper implants with the TempBase (standard procedure).

- Dimple to facilitate correct positioning
- Applicable for all implant diameters
- Available in short/long
- Multiple use
- Instruments needed: contra-angle with HXSS connection or torque wrench with surgical driver handle

Implant Driver TempBase O-Ring (spare part)

EPDM¹³, non-sterile

The Implant Driver TempBase O-Ring is used as spare part for Implant Driver TempBase.

- Multiple use

Implant Driver EV

Stainless steel¹³, PEEK⁶, non-sterile

The Implant Driver EV is used for installing the implant without using the TempBase EV.

- Multiple use
- Depth marked
- Color coded
- Instruments needed: contra-angle with HXSS connection or torque wrench with surgical driver handle

The Implant Driver EV is not the standard driver for OmniTaper implants.


Planning aids

OmniTaper Direction Indicator

Titanium¹², non-sterile

The OmniTaper Direction Indicator is used for visualizing the position and direction of the prepared osteotomy and for measuring the space between osteotomies.

- Multiple use

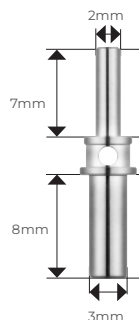

OmniTaper Direction Indicator




Order No.	68015097
-----------	----------

OmniTaper Try-In Implants

Titanium alloy², non-sterile

The OmniTaper Try-In Implant is used for intra-operative check of implant position, implant axis, congruency between implant and bone crest, and peri-implant soft tissue level.

- Design matches the OmniTaper implant
- Measure the thickness of the soft tissue
- Applicable for all implant lengths
- Multiple use
- Color coded

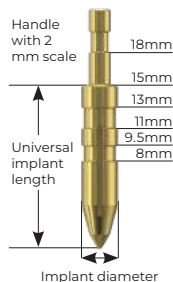

OmniTaper Try-In Implants



Ø mm	3.0	3.4	3.8	4.5	5.5
Color coding					
Order No.	68015237	68015238	68015239	68015240	68015241

Instrument Extender

Stainless steel³, non-sterile

The Instrument Extender is used for extending the length of a drill, a tap or the Implant Driver.

- Multiple use



Instrument Extender

Order No.	68015197
-----------	----------

OmniTaper Depth Gauge

Titanium alloy², non-sterile

The OmniTaper Depth Gauge is used for measuring the implant site before implant installation.

- Multiple use
- One end corresponds to the implant lengths
- One end is a measuring probe with 0-15 mm markings



OmniTaper Depth Gauge


Order No.	68015356
-----------	----------

Bone Reamer EV

Ø mm	4.0	4.6	5.2	5.8	6.4	7.0
Total length mm	26	26	26	26	26	26
Order No.	26206	26207	26208	26209	26210	26211

Bone Reamer EVStainless steel¹³, non-sterile


- The bone reamers and bone reamer guides are available from 3.0 to 5.4 and cover the use for most abutment diameters and designs.
- Laser-etched depth indication lines
- Marked with diameter
- Used for removing excess crestal bone when needed for proper abutment seating

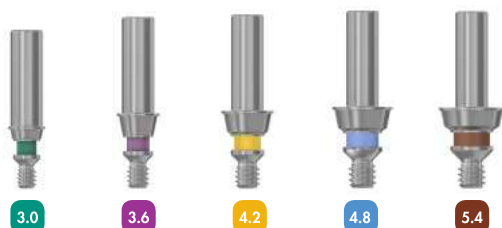Bone Reamer Guide EV

Order No.	26212	26213	26214	26215	26216
-----------	-------	-------	-------	-------	-------

Bone Reamer Guide EVStainless steel¹³, non-sterile

- Used to guide the Bone Reamer EV

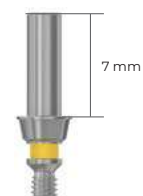**SmartFix® Guide**

Order No.	68015235
-----------	----------

SmartFix® GuideStainless steel¹³, non-sterile

- Three-piece
- Used for guiding the drilling with correct angulation


OmniTaper Radiographic Template

PC¹⁰, non-sterile

The OmniTaper Radiographic Template is a planning tool used prior to implant surgery.

- Used together with the radiographic image of the patient's jaw
- Magnification range from 1.0 to 1.5 presented in two separate sheets

OmniTaper Radiographic Template

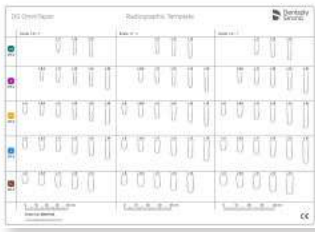

Order No.	68015132
-----------	----------

Torque Wrench

Torque Wrench EV

Stainless steel¹³, non-sterile

The Torque Wrench is used together with the Surgical Driver Handle or the Restorative Driver Handle for tightening and measuring the torque.

- Used together with the Surgical Driver Handle for implant installation and adjustment
- Used together with the Restorative Driver Handle for tightening abutments and/or bridge screws
- Multiple use

Torque Wrench EV



Order No.	25774
-----------	-------

Driver

**Torque Wrench EV
Surgical Driver Handle**

Stainless steel¹³, non-sterile

The Surgical Driver Handle is used together with the Torque Wrench.

- Multiple use

**Torque Wrench EV
Surgical Driver Handle**



Order No.	25775
-----------	-------

**Torque Wrench EV
Restorative Driver Handle**

Stainless steel¹³, non-sterile

The Restorative Driver Handle is used together with the Torque Wrench.

- Multiple use

**Torque Wrench EV
Restorative Driver Handle**



Order No.	25776	25777
-----------	-------	-------

Hex Driver manual



Total length mm	20	31
Order No.	25771	25772

Hex Driver manual

Stainless steel¹⁵, non-sterile

The Hex Driver manual is used for tightening screws, and surgical and restorative components.

Hex Driver



Total length mm	20	24	35
Order No.	68015207	68015205	68015206

Hex Driver

Stainless steel¹⁵, non-sterile

The Hex Driver is used for tightening screws, and surgical and restorative components.

OmniTaper Surgical Tray



Size mm	Length 138	Width 190	Height 61.5
Order No.	68015282		

Surgical Trays

OmniTaper Surgical Tray


PPSU⁷, silicone⁸ holders, non-sterile

The OmniTaper Surgical Tray is used during surgery and storage of OmniTaper surgical instruments.


- Instruments arranged in order of usage
- Two layers for instruments
- Grommet-less solution for easier cleaning
- Colors and icons supporting surgical procedure
- Scannable QR code for further information
- Contains no instruments

Please contact customer service to order OmniTaper Surgical Kit (surgical tray and necessary instruments)

Upper tray

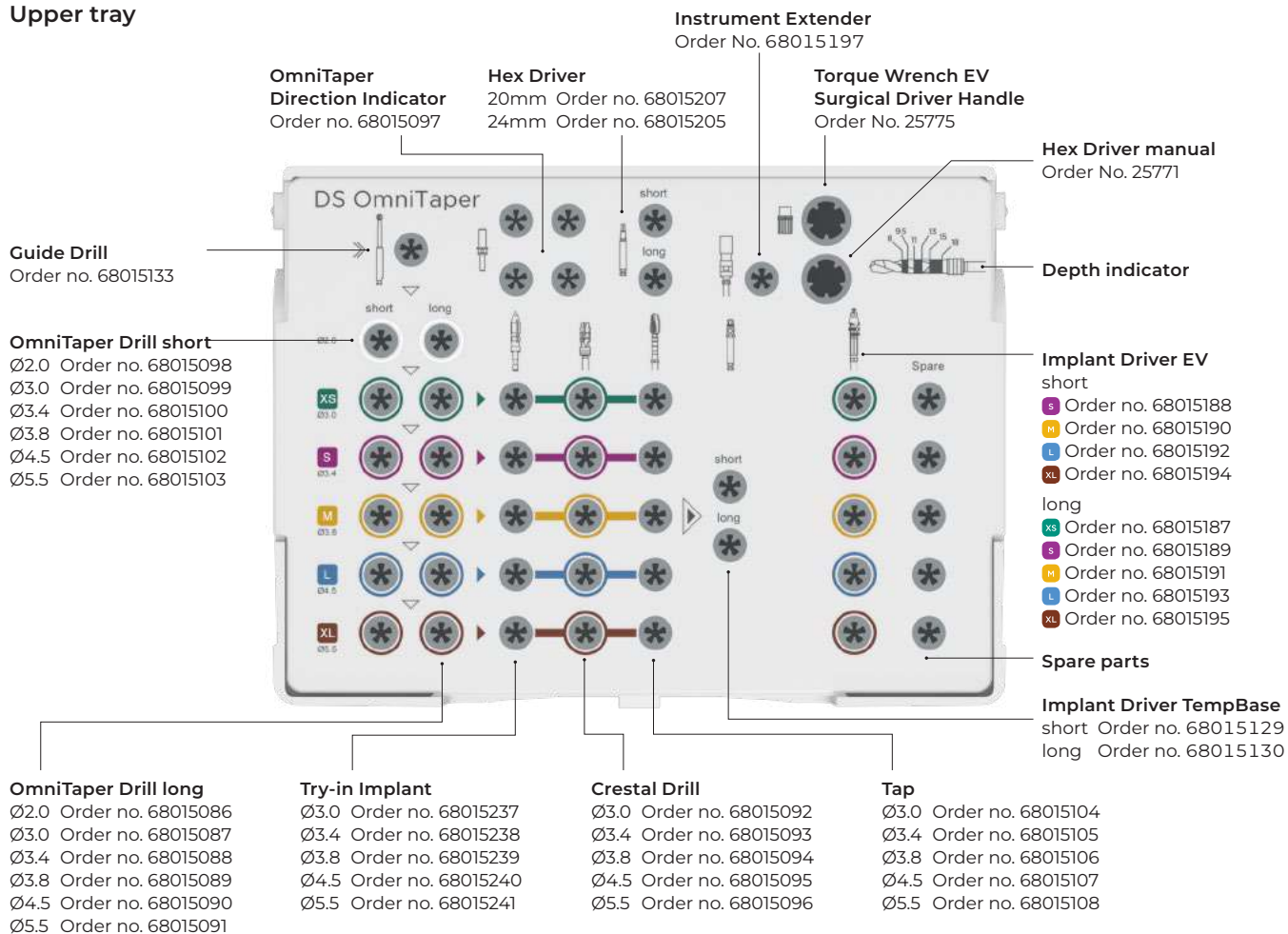

Lower tray

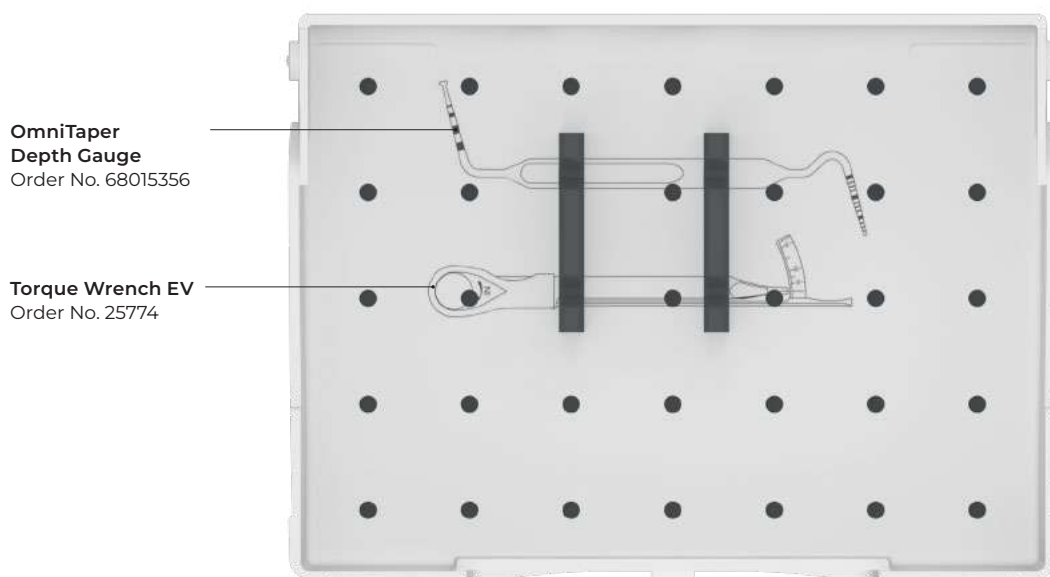

Washtray & Overlay


OmniTaper Washtray

Stainless steel³, PEEK⁴, silicon⁸, aluminum alloy⁵, non-sterile

The OmniTaper Washtray is used during surgery, for storage, machine cleaning, disinfection and subsequent steam sterilization of OmniTaper surgical instruments.

- Instrument fixation during vertical placement in washer-disinfector
- Colors and icons supporting surgical procedure
- Scannable QR code for further information
- Contains: OmniTaper Washtray and removable OmniTaper Washtray Overlay

OmniTaper Washtray

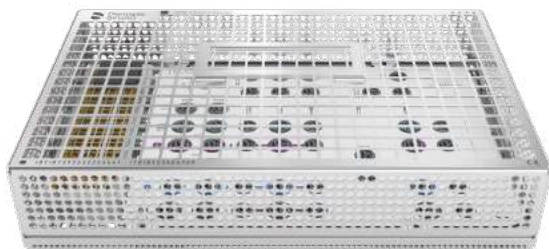

Size mm	Length 276	Width 176	Height 59
Order No.	68015271		

OmniTaper Washtray Overlay


Aluminum alloy⁶, non-sterile

The OmniTaper Washtray Overlay is used as a spare part for the OmniTaper Washtray.

- Colors and icons supporting surgical procedure

The overlay must be separated from the washtray before cleaning and re-assembled prior to sterilization. For further information, please refer to the Cleaning and sterilization instruction manual.

OmniTaper Washtray Overlay

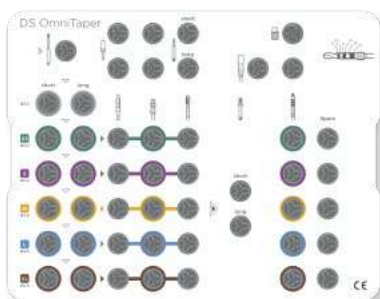

Order No.	68015272
-----------	----------

Sterilization Container & Reusable Filter

Washtray Sterilization Container

Aluminum alloy[®], Stainless steel[®], PEEK[®], Silicone[®], non-sterile

The Washtray Sterilization container is used for sterilization and storage of washtrays.

The Reusable Filter needs to be purchased separately

Washtray Sterilization Container



Order No.	68015266
-----------	----------

Spare parts	Order No.
Instrument Holder for ISO shank, PEEK-plastic with metal spring (5 pcs.)	3107 1004
Instrument Holder medium (5 pcs.)	3107 1005
Instrument Holder for Ratchet Insert, Silicone (5 pcs.)	3107 1006

Reusable Filter for Washtray Sterilization Container

PTFE[®], non-sterile

The Reusable Filter for Washtray Sterilization container is used as a spare part for Washtray Sterilization container.

- Multiple use
- Contains: 2 Reusable Filters for Washtray Sterilization Container

Reusable Filter for Washtray Sterilization Container



Order No.	68015267
-----------	----------

7. Guided Surgery Instruments

Specific instruments are available for guided implant installation of OmniTaper implants. These are marked "GS" (Guided Surgery) and can only be used together with the Simplant SAFE Guide and CEREC Guide 3. When using Simplant software or Simplant planning service, case-specific guided surgery instruments can be ordered with the Simplant SAFE Guide.

OmniTaper Mucosal Punch GS



Ø mm	3.0	3.4	3.8	4.5
Color coding				
Order No.	68017018	68017019	68017020	68017021

Punch

OmniTaper Mucosal Punch GS

Stainless steel¹³, sterile

The OmniTaper Mucosal Punch GS is used for minimally-invasive exposure of the surgical site.

- Single use
- Marked with diameter
- Color coded

OmniTaper Initial Drill GS



Ø mm	3.0	3.4	3.8	4.5
Color coding				
Order No.	68017111	68017119	68017127	68017135

Drills

OmniTaper Initial Drill GS

Stainless steel¹³, sterile

The OmniTaper Initial Drill GS is used for creating a starting point for the following drills and removal of soft tissue residue after use of the mucosa punch during guided surgery.

- Multiple use
- Marked with diameter
- Color coded

OmniTaper Drill GS



Ø mm	2.0	3.0	3.4	3.8	4.5
Color coding					
Length mm	8	8	8	8	8
Order No.	68017106	68017114	68017122	68017130	68017137

OmniTaper Drill GS

Stainless steel¹³, sterile

The OmniTaper Drill GS is used for implant site preparation up to the planned implant diameter.

- Multiple use
- Guided with the sleeve-on-drill system
- Marked with diameter
- Color coded

Drill Sleeve is needed and must be ordered separately

Drills, Crestal Drill and Tap for guided surgery can only be used with the contra-angle handpiece.

OmniTaper Drill GS

					
Ø mm	2.0	3.0	3.4	3.8	4.5
Color coding					
Length mm	9.5	9.5	9.5	9.5	9.5
Order No.	68017110	68017118	68017126	68017134	68017141

					
Ø mm	2.0	3.0	3.4	3.8	4.5
Color coding					
Length mm	11	11	11	11	11
Order No.	68017107	68017115	68017123	68017131	68017138

					
Ø mm	2.0	3.0	3.4	3.8	4.5
Color coding					
Length mm	13	13	13	13	13
Order No.	68017108	68017116	68017124	68017132	68017139



OmniTaper Drill GS


Ø mm	2.0	3.0	3.4	3.8	4.5
Color coding					
Length mm	15	15	15	15	15
Order No.	68017109	68017117	68017125	68017133	68017140

OmniTaper Crestal Drill GS

Stainless steel¹³, sterile

The OmniTaper Crestal Drill GS is used for crestal implant site preparation according to the bone quality once the planned implant diameter has been reached.

- Multiple use
- Marked with diameter
- Depth marked
- Color coded

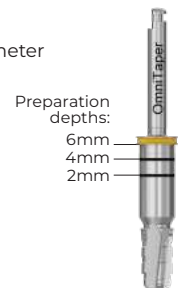

Ø mm	3.0	3.4	3.8	4.5
Color coding				
Order No.	68017100	68017101	68017102	68017103


Taps

OmniTaper Tap GS

Stainless steel¹³, sterile

The OmniTaper Tap GS is used for pre-tapping of the implant threads in cortical bone (D I).

- Multiple use
- Marked with diameter
- Depth marked (no mechanical depth stop)
- Color coded

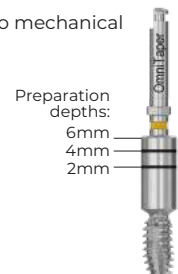

Ø mm	3.0	3.4	3.8	4.5
Color coding				
Order No.	68017014	68017015	68017016	68017017

Sleeves

OmniTaper Drill Sleeve GS

Stainless steel¹³, sterile

OmniTaper Drill Sleeve GS is used for safe guidance of the drills in the surgical guide.

- Single use
- Marked with ND (Narrow Diameter) for Ø3.0 – Ø3.8 implants or WD (Wide Diameter) for Ø4.5 implants

Must be used with the OmniTaper Drill GS

OmniTaper Drill Sleeve GS narrow diameter



Ø mm	2.0	3.0	3.4	3.8
Order No.	68017105	68017113	68017121	68017129

OmniTaper Drill Sleeve GS wide diameter



Ø mm	2.0	3.0	3.4	3.8	4.5
Order No.	68017104	68017112	68017120	68017128	68017136

Driver

OmniTaper Implant Driver TempBase GS

Stainless steel¹³, non-sterile

The OmniTaper Driver TempBase GS is used for installing OmniTaper implants with the TempBase using a contra-angle handpiece or the Torque Wrench.

- Markings on the shaft to align the abutment position and to monitor the rotational speed
- Marked with ND (Narrow Diameter) for implant Ø3.0 – Ø3.8 implants
- Marked with WD (Wide Diameter) for implant Ø4.5mm
- Sleeve is pre-mounted and can be replaced
- Torque: ≤ 50 Ncm
- Multiple use

OmniTaper Implant Driver TempBase GS short



short	ND	WD
Order No.	68017034	68017035

OmniTaper Implant Driver TempBase GS long



long	ND	WD
Order No.	68017036	68017037

Sleeves for Driver

OmniTaper Sleeve for Implant Driver TempBase GS

Stainless steel¹³, non-sterile

The OmniTaper Sleeve for Implant Driver TempBase GS is used as a spare part for the OmniTaper Implant Driver TempBase.

- Multiple use
- Marked with ND (Narrow Diameter) or WD (Wide Diameter)

OmniTaper Sleeve for Implant Driver TempBase GS



Ø	ND	WD
Order No.	68017039	68017038

Stabilization Abutment

OmniTaper Stabilization Abutment GS

Titanium alloy², non-sterile

The stabilization abutment GS for TempBase prevents lateral movement of the surgical guide between the preparation of different implant sites and the second stabilization abutment GS. It also prevents the distortion of the guide.

- Multiple use
- Narrow diameter (ND) for Ø3.0 – Ø3.8 implants
- Wide diameter (WD) for Ø4.5 implants

Positioning Aid

OmniTaper Positioning Aid EV GS

Titanium alloy², non-sterile

The OmniTaper Positioning Aid EV GS is used to secure implant analogs in the surgical guide in order to transfer the implant position, planned virtually, to a laboratory model (e.g. master cast, plaster cast).

- Multiple use

Surgical Trays

OmniTaper Surgical Tray GS

PPSU⁷, silicone⁸ holders, non-sterile

The OmniTaper Surgical Tray GS is used during surgery and storage of OmniTaper surgical instruments.

- Two layers for instruments
- Grommet-less solution for easier cleaning
- Colors and icons supporting surgical procedure
- Scannable QR code for further information
- Contains no instruments

Please contact customer service to order OmniTaper Surgical Kit GS (surgical tray and necessary instruments)

OmniTaper Stabilization Abutment GS



Ø	ND	WD
Order No.	68017043	68017045

OmniTaper Positioning Aid EV GS



Ø mm	3.0	3.4	3.8	4.5
Color coding				
Order No.	68017040	68017041	68017042	68017044

OmniTaper Surgical Tray GS



Size mm	Length 138	Width 190	Height 61.5
Order No.	68017203		

Upper tray

OmniTaper Drill GS

Ø2.0 x 8mm Order no. 68017106
 Ø2.0 x 9.5mm Order no. 68017110
 Ø2.0 x 11mm Order no. 68017107
 Ø2.0 x 13mm Order no. 68017108
 Ø2.0 x 15mm Order no. 68017109

Ø3.0 x 8mm Order no. 68017114
 Ø3.0 x 9.5mm Order no. 68017118
 Ø3.0 x 11mm Order no. 68017115
 Ø3.0 x 13mm Order no. 68017116
 Ø3.0 x 15mm Order no. 68017117

Ø3.4 x 8mm Order no. 68017122
 Ø3.4 x 9.5mm Order no. 68017126
 Ø3.4 x 11mm Order no. 68017123
 Ø3.4 x 13mm Order no. 68017124
 Ø3.4 x 15mm Order no. 68017125

Ø3.8 x 8mm Order no. 68017130
 Ø3.8 x 9.5mm Order no. 68017134
 Ø3.8 x 11mm Order no. 68017131
 Ø3.8 x 13mm Order no. 68017132
 Ø3.8 x 15mm Order no. 68017133

Ø4.5 x 8mm Order no. 68017137
 Ø4.5 x 9.5mm Order no. 68017141
 Ø4.5 x 11mm Order no. 68017138
 Ø4.5 x 13mm Order no. 68017139
 Ø4.5 x 15mm Order no. 68017140

OmniTaper Initial Drill GS


Ø3.0 Order no. 68017111
 Ø3.4 Order no. 68017119
 Ø3.8 Order no. 68017127
 Ø4.5 Order no. 68017135

OmniTaper Mucosal Punch GS

Ø3.0 Order no. 68017018
 Ø3.4 Order no. 68017019
 Ø3.8 Order no. 68017020
 Ø4.5 Order no. 68017021

Torque Wrench EV Surgical Driver Handle

Order no. 25775

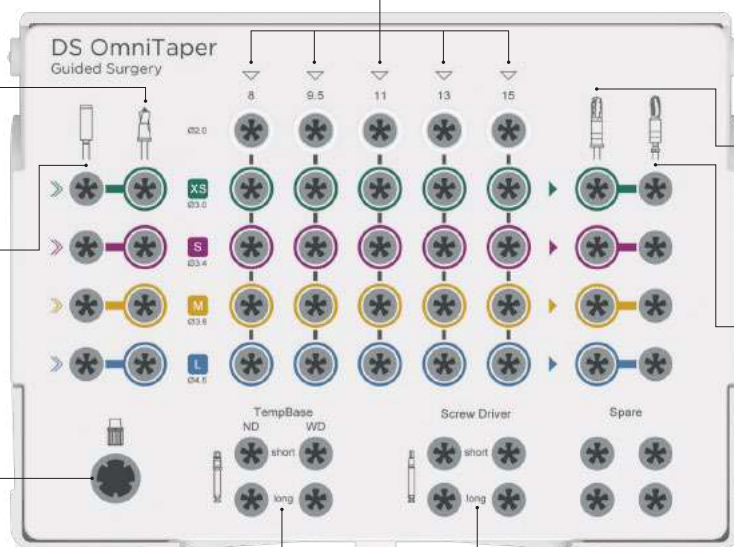

OmniTaper Crestal Drill GS

Ø3.0 Order no. 68017100
 Ø3.4 Order no. 68017101
 Ø3.8 Order no. 68017102
 Ø4.5 Order no. 68017103

OmniTaper Tap GS

Ø3.0 Order no. 68017014
 Ø3.4 Order no. 68017015
 Ø3.8 Order no. 68017016
 Ø4.5 Order no. 68017017

OmniTaper Implant Driver TempBase GS

ND short Order no. 68017034
 WD short Order no. 68017035
 ND long Order no. 68017036
 WD long Order no. 68017037

Hex Driver

20mm Order no. 68015207
 24mm Order no. 68015205


Lower tray

OmniTaper Depth Gauge


Order No. 68015356

Torque Wrench EV

Order No. 25774

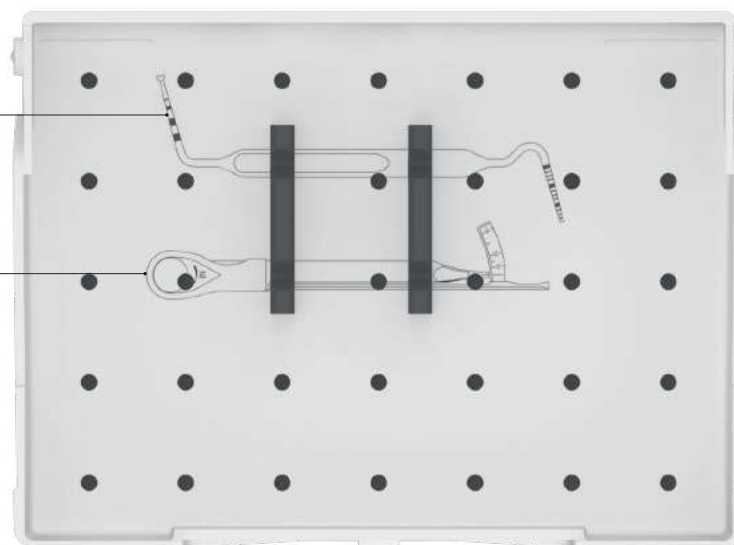

OmniTaper Washtray GS

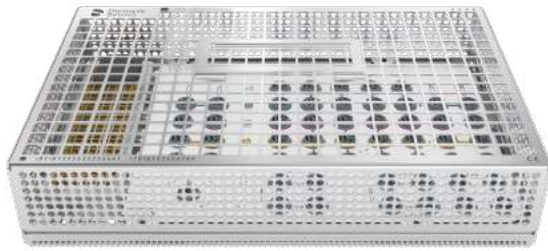

Size mm	Length 276	Width 176	Height 59
Order No.	68015354		

Spare parts

	Order No.
Instrument Holder for ISO shank, PEEK-plastic with metal spring (5 pcs.)	3107 1004
Instrument Holder medium (5 pcs.)	3107 1005
Instrument Holder for Ratchet Insert, Silicone (5 pcs.)	3107 1006


OmniTaper Washtray GS

Stainless steel¹, PEEK⁴, silicon⁸, aluminum alloy⁵, non-sterile

The OmniTaper Washtray GS is used during surgery and for storage, machine cleaning, disinfection and subsequent steam sterilization of OmniTaper surgical instruments.

- Instrument fixation during vertical placement in washer-disinfector
- Colors and icons supporting surgical procedure
- Scannable QR code for further information
- Contains: OmniTaper Washtray GS and removable OmniTaper Washtray GS Overlay

OmniTaper Washtray GS Overlay

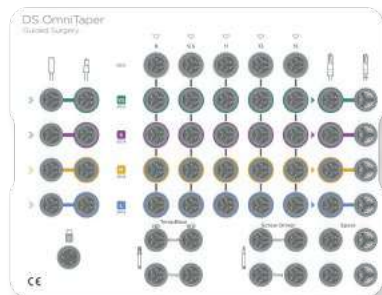

Order No.	68015355
-----------	----------

OmniTaper Washtray GS Overlay

Aluminum alloy⁵, non-sterile

The OmniTaper Washtray GS Overlay is used as a spare part for the OmniTaper Washtray GS.

- Colors and icons supporting surgical procedure

The overlay must be separated from the washtray before cleaning and re-assembled prior to sterilization. For further information, please refer to the Cleaning and sterilization instruction manual.

8. Material specifications

Index * will be found at material description for specific products.

Metals

Type	Index	Composition/Designation	Standard
Titanium	1	Ti grade 4	ASTM F67
Titanium alloy	2	Ti6Al4V ELI	ASTM F136
Stainless steel	3	Surgical stainless steel	ASTM F899
TiN Coating	4	Titanium Nitride	
Aluminum alloy	5	AlMg1	
Titanium	12	Ti grade 2	ASTM F67

Plastics

Type	Index	Composition/Designation	Standard
PEEK	6	Polyetheretherketone	
PPSU	7	Polyphenylsulfone	
Silicone	8	Polysiloxane	
ABS	9	Acrylonitrile butadiene styrene	
PC	10	Polycarbonate	
PP	11	Polypropylene	
EPDM	13	Ethylene propylene diene rubber	
PTFE	14	Polytetrafluoroethylene	

9. Index

By order number

Order No.	page	Order No.	page
25771	Hex Driver 20mm manual.....15	68013011	HealDesign EV (S) Ø5.0 4.5mm.....7
25772	Hex Driver 31mm manual.....15	68013012	HealDesign EV (S) Ø5.0 6.5mm.....7
25774	Torque Wrench EV.....14	68013013	HealDesign EV (S) Ø6.0 3.5mm.....7
25775	Torque Wrench EV Surgical Driver Handle.....14	68013014	HealDesign EV (S) Ø6.0 4.5mm.....7
25776	Torque Wrench EV Rest. Driver Handle.....14	68013015	HealDesign EV (M) Ø5.0 2.5mm.....7
25777	Torque Wrench EV Rest. Driver Handle Low.....14	68013016	HealDesign EV (M) Ø5.0 3.5mm.....7
26206	Bone Reamer EV Ø 4.0.....13	68013017	HealDesign EV (M) Ø5.0 4.5mm.....7
26207	Bone Reamer EV Ø 4.6.....13	68013018	HealDesign EV (M) Ø5.0 6.5mm.....7
26208	Bone Reamer EV Ø 5.2.....13	68013019	HealDesign EV (M) Ø6.5 2.5mm.....7
26209	Bone Reamer EV Ø 5.8.....13	68013020	HealDesign EV (M) Ø6.5 3.5mm.....7
26210	Bone Reamer EV Ø 6.4.....13	68013021	HealDesign EV (M) Ø6.5 4.5mm.....7
26211	Bone Reamer EV Ø 7.0.....13	68013022	HealDesign EV (M) Ø6.5 6.5mm.....7
26212	Bone Reamer Guide EV 3.0.....13	68013023	HealDesign EV (L) Ø5.0 2.5mm.....7
26213	Bone Reamer Guide EV 3.6.....13	68013024	HealDesign EV (L) Ø5.0 3.5mm.....7
26214	Bone Reamer Guide EV 4.2.....13	68013025	HealDesign EV (L) Ø5.0 4.5mm.....7
26215	Bone Reamer Guide EV 4.8.....13	68013026	HealDesign EV (L) Ø5.0 6.5mm.....7
26216	Bone Reamer Guide EV 5.4.....13	68013027	HealDesign EV (L) Ø6.5 2.5mm.....7
3107 1004	Instrument Holder for ISO shank.....18, 25	68013028	HealDesign EV (L) Ø6.5 3.5mm.....7
3107 1005	Instrument Holder medium.....18, 25	68013029	HealDesign EV (L) Ø6.5 4.5mm.....7
3107 1006	Instrument Holder for Ratchet Insert.....18, 25	68013030	HealDesign EV (L) Ø6.5 6.5mm.....7
68010001	Cover Screw EV (XS).....6	68013031	HealDesign EV (L) Ø7.5 3.5mm.....7
68010002	Cover Screw EV (S).....6	68013032	HealDesign EV (L) Ø7.5 4.5mm.....7
68010005	Cover Screw EV (M).....6	68013033	HealDesign EV (XL) Ø6.5 2.5mm.....7
68010006	Cover Screw EV (L).....6	68013034	HealDesign EV (XL) Ø6.5 3.5mm.....7
68010007	Cover Screw EV (XL).....6	68013035	HealDesign EV (XL) Ø6.5 4.5mm.....7
68010013	TempBase Cap (XS)-(M).....6	68013036	HealDesign EV (XL) Ø6.5 6.5mm.....7
68010014	TempBase Cap (L).....6	68013037	HealDesign EV (XL) Ø7.5 3.5mm.....7
68010015	TempBase Cap (XL).....6	68013038	HealDesign EV (XL) Ø7.5 4.5mm.....7
68011140	OmniTaper EV Ø3.0 x 11mm.....5	68013039	HealDesign EV (S) Ø5.0 3.5mm Triangular.....7
68011141	OmniTaper EV Ø3.0 x 13mm.....5	68013040	HealDesign EV (M) Ø5.0 3.5mm Triangular.....7
68011142	OmniTaper EV Ø3.0 x 15mm.....5	68013041	HealDesign EV (M) Ø6.5 3.5mm Triangular.....7
68011143	OmniTaper EV Ø3.4 x 9.5mm.....5	68013042	HealDesign EV (L) Ø5.0 3.5mm Triangular.....7
68011144	OmniTaper EV Ø3.4 x 11mm.....5	68013043	HealDesign EV (L) Ø6.5 3.5mm Triangular.....7
68011145	OmniTaper EV Ø3.4 x 13mm.....5	68013044	Healing Abutment EV (XS) 2mm.....8
68011146	OmniTaper EV Ø3.4 x 15mm.....5	68013045	Healing Abutment EV (XS) 3mm.....8
68011147	OmniTaper EV Ø3.4 x 18mm.....5	68013046	Healing Abutment EV (XS) 4mm.....8
68011148	OmniTaper EV Ø3.8 x 8mm.....5	68013047	Healing Abutment EV (XS) 6mm.....8
68011149	OmniTaper EV Ø3.8 x 9.5mm.....5	68013048	Healing Abutment EV (S) 2mm.....8
68011150	OmniTaper EV Ø3.8 x 11mm.....5	68013049	Healing Abutment EV (S) 3mm.....8
68011151	OmniTaper EV Ø3.8 x 13mm.....5	68013050	Healing Abutment EV (S) 4mm.....8
68011152	OmniTaper EV Ø3.8 x 15mm.....5	68013051	Healing Abutment EV (S) 6mm.....8
68011153	OmniTaper EV Ø3.8 x 18mm.....5	68013052	Healing Abutment EV (M) 2mm.....8
68011154	OmniTaper EV Ø4.5 x 8mm.....5	68013053	Healing Abutment EV (M) 3mm.....8
68011155	OmniTaper EV Ø4.5 x 9.5mm.....5	68013054	Healing Abutment EV (M) 4mm.....8
68011156	OmniTaper EV Ø4.5 x 11mm.....5	68013055	Healing Abutment EV (M) 6mm.....8
68011157	OmniTaper EV Ø4.5 x 13mm.....5	68013056	Healing Abutment EV (L) 2mm.....8
68011158	OmniTaper EV Ø4.5 x 15mm.....5	68013057	Healing Abutment EV (L) 3mm.....8
68011159	OmniTaper EV Ø4.5 x 18mm.....5	68013058	Healing Abutment EV (L) 4mm.....8
68011160	OmniTaper EV Ø5.5 x 8mm.....5	68013059	Healing Abutment EV (L) 6mm.....8
68011161	OmniTaper EV Ø5.5 x 9.5mm.....5	68013060	Healing Abutment EV (XL) 2mm.....8
68011162	OmniTaper EV Ø5.5 x 11mm.....5	68013061	Healing Abutment EV (XL) 3mm.....8
68011163	OmniTaper EV Ø5.5 x 13mm.....5	68013062	Healing Abutment EV (XL) 4mm.....8
68011164	OmniTaper EV Ø5.5 x 15mm.....5	68013063	Healing Abutment EV (XL) 6mm.....8
68013001	HealDesign EV (XS) Ø3.5 2.5mm.....6	68015080	OmniTaper Drill SP Ø2.0 short.....10
68013002	HealDesign EV (XS) Ø3.5 3.5mm.....6	68015081	OmniTaper Drill SP Ø3.0 short.....10
68013003	HealDesign EV (XS) Ø3.5 4.5mm.....6	68015082	OmniTaper Drill SP Ø3.4 short.....10
68013004	HealDesign EV (XS) Ø3.5 6.5mm.....6	68015083	OmniTaper Drill SP Ø3.8 short.....10
68013005	HealDesign EV (S) Ø4.0 2.5mm.....7	68015084	OmniTaper Drill SP Ø4.5 short.....10
68013006	HealDesign EV (S) Ø4.0 3.5mm.....7	68015085	OmniTaper Drill SP Ø5.5 short.....10
68013007	HealDesign EV (S) Ø4.0 4.5mm.....7	68015086	OmniTaper Drill Ø2.0 long.....9
68013008	HealDesign EV (S) Ø4.0 6.5mm.....7	68015087	OmniTaper Drill Ø3.0 long.....9
68013009	HealDesign EV (S) Ø5.0 2.5mm.....7	68015088	OmniTaper Drill Ø3.4 long.....9
68013010	HealDesign EV (S) Ø5.0 3.5mm.....7	68015089	OmniTaper Drill Ø3.8 long.....9

Order No.	page
68015090	OmniTaper Drill Ø4.5 long.....9
68015091	OmniTaper Drill Ø5.5 long.....9
68015092	OmniTaper Crestal Drill Ø3.0.....10
68015093	OmniTaper Crestal Drill Ø3.4.....10
68015094	OmniTaper Crestal Drill Ø3.8.....10
68015095	OmniTaper Crestal Drill Ø4.5.....10
68015096	OmniTaper Crestal Drill Ø5.5.....10
68015097	OmniTaper Direction Indicator.....12
68015098	OmniTaper Drill Ø2.0 short.....9
68015099	OmniTaper Drill Ø3.0 short.....9
68015100	OmniTaper Drill Ø3.4 short.....9
68015101	OmniTaper Drill Ø3.8 short.....9
68015102	OmniTaper Drill Ø4.5 short.....9
68015103	OmniTaper Drill Ø5.5 short.....9
68015104	OmniTaper Tap Ø3.0.....10
68015105	OmniTaper Tap Ø3.4.....10
68015106	OmniTaper Tap Ø3.8.....10
68015107	OmniTaper Tap Ø4.5.....10
68015108	OmniTaper Tap Ø5.5.....10
68015129	Implant Driver TempBase short.....11
68015130	Implant Driver TempBase long.....11
68015131	Implant Driver TempBase O-Ring.....11
68015132	OmniTaper Radiographic Template.....14
68015133	Guide Drill.....9
68015134	Precision Drill.....9
68015187	Implant Driver EV (XS) long.....11
68015188	Implant Driver EV (S) short.....11
68015189	Implant Driver EV (S) long.....11
68015190	Implant Driver EV (M) short.....11
68015191	Implant Driver EV (M) long.....11
68015192	Implant Driver EV (L) short.....11
68015193	Implant Driver EV (L) long.....11
68015194	Implant Driver EV (XL) short.....11
68015195	Implant Driver EV (XL) long.....11
68015197	Instrument Extender.....12
68015205	Hex Driver 24mm.....15
68015206	Hex Driver 35mm.....15
68015207	Hex Driver 20mm.....15
68015235	SmartFix Guide.....13
68015237	OmniTaper Try-in Implant Ø3.0.....12
68015238	OmniTaper Try-in Implant Ø3.4.....12
68015239	OmniTaper Try-in Implant Ø3.8.....12
68015240	OmniTaper Try-in Implant Ø4.5.....12
68015241	OmniTaper Try-in Implant Ø5.5.....12
68015266	Washtray Sterilization Container.....18
68015267	Washtray Sterilization Container Filter 2x.....18
68015271	OmniTaper Washtray.....17
68015272	OmniTaper Washtray Overlay.....17
68015282	OmniTaper Surgical Tray.....15
68015354	OmniTaper Washtray GS.....25
68015355	OmniTaper Washtray Overlay GS.....25
68015356	OmniTaper Depth Gauge.....12
68017014	OmniTaper Tap Ø3.0 GS.....21
68017015	OmniTaper Tap Ø3.4 GS.....21
68017016	OmniTaper Tap Ø3.8 GS.....21
68017017	OmniTaper Tap Ø4.5 GS.....21
68017018	OmniTaper Mucosal Punch Ø3.0 GS.....19
68017019	OmniTaper Mucosal Punch Ø3.4 GS.....19

Order No.	page
68017020	OmniTaper Mucosal Punch Ø3.8 GS.....19
68017021	OmniTaper Mucosal Punch Ø4.5 GS.....19
68017034	OmniTaper Implant Driver TempBase ND short GS.....22
68017035	OmniTaper Implant Driver TempBase WD short GS.....22
68017036	OmniTaper Implant Driver TempBase ND long GS.....22
68017037	OmniTaper Implant Driver TempBase WD long GS.....22
68017038	OmniTaper Sleeve for Implant Driver TempBase WD GS.....22
68017039	OmniTaper Sleeve for Implant Driver TempBase ND GS.....22
68017040	OmniTaper Positioning Aid EV (XS) GS.....23
68017041	OmniTaper Positioning Aid EV (S) GS.....23
68017042	OmniTaper Positioning Aid EV (M) GS.....23
68017043	OmniTaper Stabilization Abutment ND GS.....23
68017044	OmniTaper Positioning Aid EV (L) GS.....23
68017045	OmniTaper Stabilization Abutment WD GS.....23
68017100	OmniTaper Crestal Drill Ø3.0 GS.....21
68017101	OmniTaper Crestal Drill Ø3.4 GS.....21
68017102	OmniTaper Crestal Drill Ø3.8 GS.....21
68017103	OmniTaper Crestal Drill Ø4.5 GS.....21
68017104	OmniTaper Drill Sleeve WD for Drill Ø2.0 GS.....22
68017105	OmniTaper Drill Sleeve ND for Drill Ø2.0 GS.....22
68017106	OmniTaper Drill Ø2.0 x 8mm GS.....19
68017107	OmniTaper Drill Ø2.0 x 11mm GS.....18
68017108	OmniTaper Drill Ø2.0 x 13mm GS.....18
68017109	OmniTaper Drill Ø2.0 x 15mm GS.....21
68017110	OmniTaper Drill Ø2.0 x 9.5mm GS.....20
68017111	OmniTaper Initial Drill Ø3.0 GS.....19
68017112	OmniTaper Drill Sleeve WD for Drill Ø3.0 GS.....22
68017113	OmniTaper Drill Sleeve ND for Drill Ø3.0 GS.....22
68017114	OmniTaper Drill Ø3.0 x 8mm GS.....19
68017115	OmniTaper Drill Ø3.0 x 11mm GS.....20
68017116	OmniTaper Drill Ø3.0 x 13mm GS.....20
68017117	OmniTaper Drill Ø3.0 x 15mm GS.....21
68017118	OmniTaper Drill Ø3.0 x 9.5mm GS.....20
68017119	OmniTaper Initial Drill Ø3.4 GS.....19
68017120	OmniTaper Drill Sleeve WD for Drill Ø3.4 GS.....22
68017121	OmniTaper Drill Sleeve ND for Drill Ø3.4 GS.....22
68017122	OmniTaper Drill Ø3.4 x 8mm GS.....19
68017123	OmniTaper Drill Ø3.4 x 11mm GS.....20
68017124	OmniTaper Drill Ø3.4 x 13mm GS.....20
68017125	OmniTaper Drill Ø3.4 x 15mm GS.....21
68017126	OmniTaper Drill Ø3.4 x 9.5mm GS.....20
68017127	OmniTaper Initial Drill Ø3.8 GS.....19
68017128	OmniTaper Drill Sleeve WD for Drill Ø3.8 GS.....22
68017129	OmniTaper Drill Sleeve ND for Drill Ø3.8 GS.....22
68017130	OmniTaper Drill Ø3.8 x 8mm GS.....19
68017131	OmniTaper Drill Ø3.8 x 11mm GS.....20
68017132	OmniTaper Drill Ø3.8 x 13mm GS.....20
68017133	OmniTaper Drill Ø3.8 x 15mm GS.....21
68017134	OmniTaper Drill Ø3.8 x 9.5mm GS.....20
68017135	OmniTaper Initial Drill Ø4.5 GS.....19
68017136	OmniTaper Drill Sleeve WD for Drill Ø4.5 GS.....22
68017137	OmniTaper Drill Ø4.5 x 8mm GS.....19
68017138	OmniTaper Drill Ø4.5 x 11mm GS.....20
68017139	OmniTaper Drill Ø4.5 x 13mm GS.....20
68017140	OmniTaper Drill Ø4.5 x 15mm GS.....21
68017141	OmniTaper Drill Ø4.5 x 9.5mm GS.....20
68017203	OmniTaper Surgical Tray GS.....23

Alphabetical order

Order No.	page	Order No.	page
26206	Bone Reamer EV Ø 4.0.....13	68013059	Healing Abutment EV (L) 6mm.....8
26207	Bone Reamer EV Ø 4.6.....13	68013052	Healing Abutment EV (M) 2mm.....8
26208	Bone Reamer EV Ø 5.2.....13	68013053	Healing Abutment EV (M) 3mm.....8
26209	Bone Reamer EV Ø 5.8.....13	68013054	Healing Abutment EV (M) 4mm.....8
26210	Bone Reamer EV Ø 6.4.....13	68013055	Healing Abutment EV (M) 6mm.....8
26211	Bone Reamer EV Ø 7.0.....13	68013048	Healing Abutment EV (S) 2mm.....8
26212	Bone Reamer Guide EV 3.0.....13	68013049	Healing Abutment EV (S) 3mm.....8
26213	Bone Reamer Guide EV 3.6.....13	68013050	Healing Abutment EV (S) 4mm.....8
26214	Bone Reamer Guide EV 4.2.....13	68013051	Healing Abutment EV (S) 6mm.....8
26215	Bone Reamer Guide EV 4.8.....13	68013060	Healing Abutment EV (XL) 2mm.....8
26216	Bone Reamer Guide EV 5.4.....13	68013061	Healing Abutment EV (XL) 3mm.....8
68010006	Cover Screw EV (L).....6	68013062	Healing Abutment EV (XL) 4mm.....8
68010005	Cover Screw EV (M).....6	68013063	Healing Abutment EV (XL) 6mm.....8
68010002	Cover Screw EV (S).....6	68013044	Healing Abutment EV (XS) 2mm.....8
68010007	Cover Screw EV (XL).....6	68013045	Healing Abutment EV (XS) 3mm.....8
68010001	Cover Screw EV (XS).....6	68013046	Healing Abutment EV (XS) 4mm.....8
68015133	Guide Drill.....9	68013047	Healing Abutment EV (XS) 6mm.....8
68013023	HealDesign EV (L) Ø5.0 2.5mm.....7	68015207	Hex Driver 20mm.....15
68013024	HealDesign EV (L) Ø5.0 3.5mm.....7	25771	Hex Driver 20mm manual.....15
68013042	HealDesign EV (L) Ø5.0 3.5mm Triangular.....7	68015205	Hex Driver 24mm.....15
68013025	HealDesign EV (L) Ø5.0 4.5mm.....7	25772	Hex Driver 31mm manual.....15
68013026	HealDesign EV (L) Ø5.0 6.5mm.....7	68015206	Hex Driver 35mm.....15
68013027	HealDesign EV (L) Ø6.5 2.5mm.....7	68015193	Implant Driver EV (L) long.....11
68013028	HealDesign EV (L) Ø6.5 3.5mm.....7	68015192	Implant Driver EV (L) short.....11
68013043	HealDesign EV (L) Ø6.5 3.5mm Triangular.....7	68015191	Implant Driver EV (M) long.....11
68013029	HealDesign EV (L) Ø6.5 4.5mm.....7	68015190	Implant Driver EV (M) short.....11
68013030	HealDesign EV (L) Ø6.5 6.5mm.....7	68015189	Implant Driver EV (S) long.....11
68013031	HealDesign EV (L) Ø7.5 3.5mm.....7	68015188	Implant Driver EV (S) short.....11
68013032	HealDesign EV (L) Ø7.5 4.5mm.....7	68015195	Implant Driver EV (XL) long.....11
68013015	HealDesign EV (M) Ø5.0 2.5mm.....7	68015194	Implant Driver EV (XL) short.....11
68013016	HealDesign EV (M) Ø5.0 3.5mm.....7	68015187	Implant Driver EV (XS) long.....11
68013040	HealDesign EV (M) Ø5.0 3.5mm Triangular.....7	68015130	Implant Driver TempBase long.....11
68013017	HealDesign EV (M) Ø5.0 4.5mm.....7	68015131	Implant Driver TempBase O-Ring.....11
68013018	HealDesign EV (M) Ø5.0 6.5mm.....7	68015129	Implant Driver TempBase short.....11
68013019	HealDesign EV (M) Ø6.5 2.5mm.....7	68015197	Instrument Extender.....12
68013020	HealDesign EV (M) Ø6.5 3.5mm.....7	3107 1004	Instrument Holder for ISO shank.....18, 25
68013041	HealDesign EV (M) Ø6.5 3.5mm Triangular.....7	3107 1006	Instrument Holder for Ratchet Insert.....18, 25
68013021	HealDesign EV (M) Ø6.5 4.5mm.....7	3107 1005	Instrument Holder medium.....18, 25
68013022	HealDesign EV (M) Ø6.5 6.5mm.....7	68015092	OmniTaper Crestal Drill Ø3.0.....10
68013005	HealDesign EV (S) Ø4.0 2.5mm.....7	68017100	OmniTaper Crestal Drill Ø3.0 GS.....21
68013006	HealDesign EV (S) Ø4.0 3.5mm.....7	68015093	OmniTaper Crestal Drill Ø3.4.....10
68013007	HealDesign EV (S) Ø4.0 4.5mm.....7	68017101	OmniTaper Crestal Drill Ø3.4 GS.....21
68013008	HealDesign EV (S) Ø4.0 6.5mm.....7	68015094	OmniTaper Crestal Drill Ø3.8.....10
68013009	HealDesign EV (S) Ø5.0 2.5mm.....7	68017102	OmniTaper Crestal Drill Ø3.8 GS.....21
68013010	HealDesign EV (S) Ø5.0 3.5mm.....7	68015095	OmniTaper Crestal Drill Ø4.5.....10
68013039	HealDesign EV (S) Ø5.0 3.5mm Triangular.....7	68017103	OmniTaper Crestal Drill Ø4.5 GS.....21
68013011	HealDesign EV (S) Ø5.0 4.5mm.....7	68015096	OmniTaper Crestal Drill Ø5.5.....10
68013012	HealDesign EV (S) Ø5.0 6.5mm.....7	68015356	OmniTaper Depth Gauge.....12
68013013	HealDesign EV (S) Ø6.0 3.5mm.....7	68015097	OmniTaper Direction Indicator.....12
68013014	HealDesign EV (S) Ø6.0 4.5mm.....7	68015086	OmniTaper Drill Ø2.0 long.....9
68013033	HealDesign EV (XL) Ø6.5 2.5mm.....7	68015098	OmniTaper Drill Ø2.0 short.....9
68013034	HealDesign EV (XL) Ø6.5 3.5mm.....7	68017106	OmniTaper Drill Ø2.0 x 8mm GS.....19
68013035	HealDesign EV (XL) Ø6.5 4.5mm.....7	68017110	OmniTaper Drill Ø2.0 x 9.5mm GS.....20
68013036	HealDesign EV (XL) Ø6.5 6.5mm.....7	68017107	OmniTaper Drill Ø2.0 x 11mm GS.....20
68013037	HealDesign EV (XL) Ø7.5 3.5mm.....7	68017108	OmniTaper Drill Ø2.0 x 13mm GS.....20
68013038	HealDesign EV (XL) Ø7.5 4.5mm.....7	68017109	OmniTaper Drill Ø2.0 x 15mm GS.....21
68013001	HealDesign EV (XS) Ø3.5 2.5mm.....6	68015087	OmniTaper Drill Ø3.0 long.....9
68013002	HealDesign EV (XS) Ø3.5 3.5mm.....6	68015099	OmniTaper Drill Ø3.0 short.....9
68013003	HealDesign EV (XS) Ø3.5 4.5mm.....6	68017114	OmniTaper Drill Ø3.0 x 8mm GS.....19
68013004	HealDesign EV (XS) Ø3.5 6.5mm.....6	68017118	OmniTaper Drill Ø3.0 x 9.5mm GS.....20
68013056	Healing Abutment EV (L) 2mm.....8	68017115	OmniTaper Drill Ø3.0 x 11mm GS.....20
68013057	Healing Abutment EV (L) 3mm.....8	68017116	OmniTaper Drill Ø3.0 x 13mm GS.....20
68013058	Healing Abutment EV (L) 4mm.....8	68017117	OmniTaper Drill Ø3.0 x 15mm GS.....21

Order No.	page
68015088	OmniTaper Drill Ø3.4 long.....9
68015100	OmniTaper Drill Ø3.4 short.....9
68017122	OmniTaper Drill Ø3.4 x 8mm GS.....19
68017126	OmniTaper Drill Ø3.4 x 9.5mm GS.....20
68017123	OmniTaper Drill Ø3.4 x 11mm GS.....20
68017124	OmniTaper Drill Ø3.4 x 13mm GS.....20
68017125	OmniTaper Drill Ø3.4 x 15mm GS.....21
68015089	OmniTaper Drill Ø3.8 long.....9
68015101	OmniTaper Drill Ø3.8 short.....9
68017130	OmniTaper Drill Ø3.8 x 8mm GS.....19
68017134	OmniTaper Drill Ø3.8 x 9.5mm GS.....20
68017131	OmniTaper Drill Ø3.8 x 11mm GS.....20
68017132	OmniTaper Drill Ø3.8 x 13mm GS.....20
68017133	OmniTaper Drill Ø3.8 x 15mm GS.....21
68015090	OmniTaper Drill Ø4.5 long.....9
68015102	OmniTaper Drill Ø4.5 short.....9
68017137	OmniTaper Drill Ø4.5 x 8mm GS.....19
68017141	OmniTaper Drill Ø4.5 x 9.5mm GS.....20
68017138	OmniTaper Drill Ø4.5 x 11mm GS.....20
68017139	OmniTaper Drill Ø4.5 x 13mm GS.....20
68017140	OmniTaper Drill Ø4.5 x 15mm GS.....21
68015091	OmniTaper Drill Ø5.5 long.....9
68015103	OmniTaper Drill Ø5.5 short.....9
68017105	OmniTaper Drill Sleeve ND for Drill Ø2.0 GS.....22
68017113	OmniTaper Drill Sleeve ND for Drill Ø3.0 GS.....22
68017121	OmniTaper Drill Sleeve ND for Drill Ø3.4 GS.....22
68017129	OmniTaper Drill Sleeve ND for Drill Ø3.8 GS.....22
68017104	OmniTaper Drill Sleeve WD for Drill Ø2.0 GS.....22
68017112	OmniTaper Drill Sleeve WD for Drill Ø3.0 GS.....22
68017120	OmniTaper Drill Sleeve WD for Drill Ø3.4 GS.....22
68017128	OmniTaper Drill Sleeve WD for Drill Ø3.8 GS.....22
68017136	OmniTaper Drill Sleeve WD for Drill Ø4.5 GS.....22
68015080	OmniTaper Drill SP Ø2.0 short.....10
68015081	OmniTaper Drill SP Ø3.0 short.....10
68015082	OmniTaper Drill SP Ø3.4 short.....10
68015083	OmniTaper Drill SP Ø3.8 short.....10
68015084	OmniTaper Drill SP Ø4.5 short.....10
68015085	OmniTaper Drill SP Ø5.5 short.....10
68011140	OmniTaper EV Ø3.0 x 11mm.....5
68011141	OmniTaper EV Ø3.0 x 13mm.....5
68011142	OmniTaper EV Ø3.0 x 15mm.....5
68011143	OmniTaper EV Ø3.4 x 9.5mm.....5
68011144	OmniTaper EV Ø3.4 x 11mm.....5
68011145	OmniTaper EV Ø3.4 x 13mm.....5
68011146	OmniTaper EV Ø3.4 x 15mm.....5
68011147	OmniTaper EV Ø3.4 x 18mm.....5
68011148	OmniTaper EV Ø3.8 x 8mm.....5
68011149	OmniTaper EV Ø3.8 x 9.5mm.....5
68011150	OmniTaper EV Ø3.8 x 11mm.....5
68011151	OmniTaper EV Ø3.8 x 13mm.....5
68011152	OmniTaper EV Ø3.8 x 15mm.....5
68011153	OmniTaper EV Ø3.8 x 18mm.....5
68011154	OmniTaper EV Ø4.5 x 8mm.....5
68011155	OmniTaper EV Ø4.5 x 9.5mm.....5
68011156	OmniTaper EV Ø4.5 x 11mm.....5
68011157	OmniTaper EV Ø4.5 x 13mm.....5
68011158	OmniTaper EV Ø4.5 x 15mm.....5
68011159	OmniTaper EV Ø4.5 x 18mm.....5

Order No.	page
68011160	OmniTaper EV Ø5.5 x 8mm.....5
68011161	OmniTaper EV Ø5.5 x 9.5mm.....5
68011162	OmniTaper EV Ø5.5 x 11mm.....5
68011163	OmniTaper EV Ø5.5 x 13mm.....5
68011164	OmniTaper EV Ø5.5 x 15mm.....5
68017036	OmniTaper Implant Driver TempBase ND long GS.....22
68017034	OmniTaper Implant Driver TempBase ND short GS.....22
68017037	OmniTaper Implant Driver TempBase WD long GS.....22
68017035	OmniTaper Implant Driver TempBase WD short GS.....22
68017111	OmniTaper Initial Drill Ø3.0 GS.....18
68017119	OmniTaper Initial Drill Ø3.4 GS.....18
68017127	OmniTaper Initial Drill Ø3.8 GS.....18
68017135	OmniTaper Initial Drill Ø4.5 GS.....18
68017018	OmniTaper Mucosal Punch Ø3.0 GS.....18
68017019	OmniTaper Mucosal Punch Ø3.4 GS.....18
68017020	OmniTaper Mucosal Punch Ø3.8 GS.....18
68017021	OmniTaper Mucosal Punch Ø4.5 GS.....18
68017044	OmniTaper Positioning Aid EV (L) GS.....23
68017042	OmniTaper Positioning Aid EV (M) GS.....23
68017041	OmniTaper Positioning Aid EV (S) GS.....23
68017040	OmniTaper Positioning Aid EV (XS) GS.....23
68015132	OmniTaper Radiographic Template.....14
68017039	OmniTaper Sleeve for Implant Driver TempBase ND GS.....22
68017038	OmniTaper Sleeve for Implant Driver TempBase WD GS.....22
68017043	OmniTaper Stabilization Abutment ND GS.....23
68017045	OmniTaper Stabilization Abutment WD GS.....23
68015282	OmniTaper Surgical Tray.....15
68017203	OmniTaper Surgical Tray GS.....23
68015104	OmniTaper Tap Ø3.0.....10
68017014	OmniTaper Tap Ø3.0 GS.....21
68015105	OmniTaper Tap Ø3.4.....10
68017015	OmniTaper Tap Ø3.4 GS.....21
68015106	OmniTaper Tap Ø3.8.....10
68017016	OmniTaper Tap Ø3.8 GS.....21
68015107	OmniTaper Tap Ø4.5.....10
68017017	OmniTaper Tap Ø4.5 GS.....21
68015108	OmniTaper Tap Ø5.5.....10
68015237	OmniTaper Try-in Implant Ø3.0.....12
68015238	OmniTaper Try-in Implant Ø3.4.....12
68015239	OmniTaper Try-in Implant Ø3.8.....12
68015240	OmniTaper Try-in Implant Ø4.5.....12
68015241	OmniTaper Try-in Implant Ø5.5.....12
68015271	OmniTaper Washtray.....17
68015354	OmniTaper Washtray GS.....25
68015272	OmniTaper Washtray Overlay.....17
68015355	OmniTaper Washtray Overlay GS.....25
68015134	Precision Drill.....9
68015235	SmartFix Guide.....13
68010014	TempBase Cap (L).....6
68010015	TempBase Cap (XL).....6
68010013	TempBase Cap (XS)-(M).....6
25774	Torque Wrench EV.....14
25776	Torque Wrench EV Rest. Driver Handle.....14
25775	Torque Wrench EV Surgical Driver Handle.....14
25777	Torque Wrench EV Rest. Driver Handle Low.....14
68015266	Washtray Sterilization Container.....18
68015267	Washtray Sterilization Container Filter 2x.....18

Notes

[illegible]

